MMT-8007 - DIGITAL ALARM-THERMOMETER WITH PROBE

SPECIFICATION

1.Measuring temperature range:
- External probe Temperature Range: -50~+70˚C/-58~+158˚F
- Internal sensor Temperature Range: 0~+50˚C/+32~+122˚F

2.Temperature accuracy: ±1℃ (±2℉)

3.Temperature display resolution: 0.1℃ (0.1℉)

4.MAX/MIN temperature memory records

5.Operating voltage:1.5V, 1 PCS AAA battery

6.Total Dimension:110mm*70mm*19.5mm

7.Net weight: pprox.82g

FEATURES

1.LCD screen shows indoor & outdoor temperature.

2.Preset alert temperature.

3.Stand for shelf or hanging hole at back.

4.Max/Min temperature memory records.

5.The MMT-8007A has 2 channels:
- The top display labelled “IN” is the internal sensor in the display
- The bottom display labelled “OUT” is the external wired sensor

BUTTONS & OPERATIONS

1. MIN/MAX Button:

a) Press MIN/MAX key once to display Max temp on IN/OUT channels
b) Press MIN/MAX key again to display Min temp on IN/OUT channels
c) Press MIN/MAX key again to return to main screen

2. °C/°F Button: Alt between C/F

3. RESET Button:

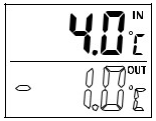
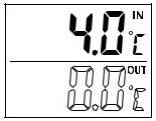
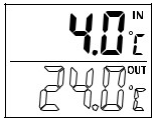
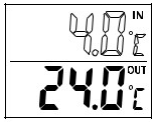
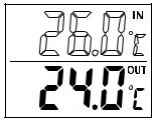
a) Press once = clears min/max temperatures
b) Press twice = reset Min/Max range to factor default (-1℃ to 4℃)

4. ALERT/SET Button:

a) Press once to turn on/off Alert
b) Hold 1 second to enter alert temperature range
c) OUT temperature flashes. Press C/F to increase OUT temperature to 2.0°C
d) Press Alert/Set again. Press C/F to increase IN temperature to 8°C
e) Press Alert/Set to set alarm range at 2° to 8°C as per instructions above (this range can be customized as per user’s need)

**Note: Make sure to set Alert ON. Place it in refrigerator to activate alarm function.*

Scigiene
1295 Morningside Avenue, Unit 16-18
Scarborough, ON M1B 4Z4 Canada
Phone: 416-261-4865 Fax: 416-261-7879
www.scigiene.com

ALERT SET	Applies to External Probe (OUT)	
ALERT SET	HOLDING 1 SECOND FOR SETTING ALERT TEMPERATURE RANGE to see default setting Press once to adjust “OUT” (MAX) limit	
C/F	PRESS ONCE ADDING ONE DEG C at a time	
C/F	HOLDING THE KEY will advance temperature rapidly	
ALERT SET	Press 2X to adjust “IN” (MIN) limit	
C/F	PRESS ONCE ADDING ONE DEG C at a time. Holding the key will advance temperature rapidly.	
ALERT SET	Once activated as per note in section 4, buzzer on 5 seconds and off 55 seconds when the Outdoor Temp. is between (+24degC~ +26degC.)	

Scigiene
1295 Morningside Avenue, Unit 16-18
Scarborough, ON M1B 4Z4 Canada
Phone: 416-261-4865 Fax: 416-261-7879
www.scigiene.com